



Unlock the Latest in Skin Tightening & Lifting

6.78 MHz Monopolar RF
Continuous Water Cooling
Hidden Edge Technology
Painless & Effective



VOLNEWMER

| 6.78MHz Monopolar Radiofrequency

Transform your skincare with Volnewmer, an FDA-approved innovation for treating skin laxity on the face and body.

With advanced, non-invasive RF technology, Volnewmer delivers visible and effective results with minimal discomfort, powered by five optimized technologies for precise energy delivery.

Delivers **TTP**

(**Targeted Thermal Precision RF**) evenly to the dermis while protecting the epidermis, ensuring safe and effective tightening and lifting results.

Ensures a comfortable treatment experience for patients without pain, with noticeable results and minimal discomfort

With flexible tips and impedance monitoring, full shots deliver immediate and lasting tightening/lifting outcomes



Easy to combine with other treatments
Effective for both face and body, suitable for all skin types.



Unlocking the Science Within


The hidden science in Volnewmer, drives superior clinical outcomes

1 Delivers full double depth shots for deep thermal transformation and collagen generation.

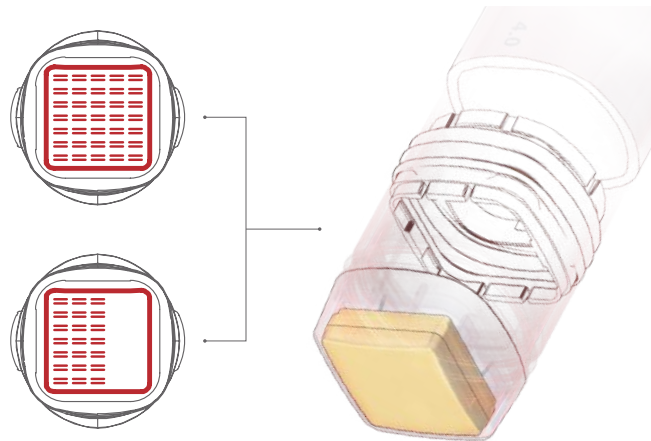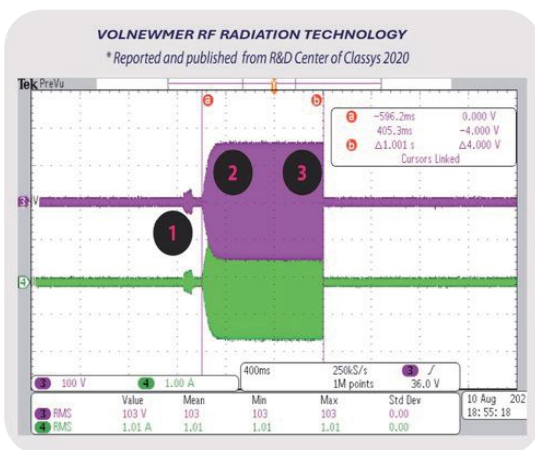

- ✓ **Long-lasting clinical outcomes**, through thermal energy delivery to the reticular dermis
- ✓ A single pulse delivers consistent RF energy without a treatment gap, maintaining long-lasting thermal energy that stimulates collagen generation in both the papillary dermis and reticular dermis

2 The continuous water-cooling system maintains safe, consistent epidermal temperatures while delivering optimal energy deep into the dermis for effective results.

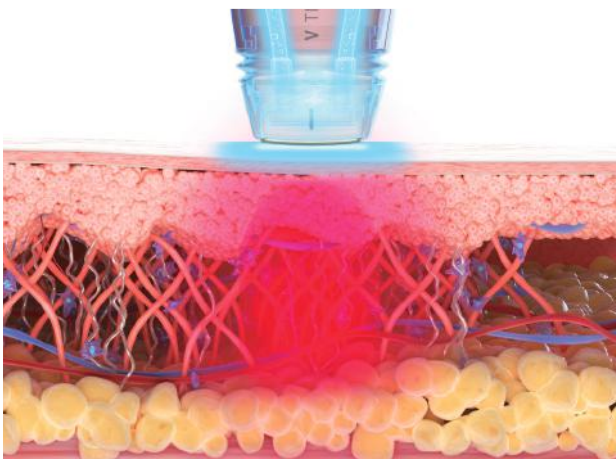

- ✓ **No Pain**
Ensures minimal to no pain during treatment
- ✓ **Deeper Energy Delivery**
Targets both papillary dermis and reticular dermis
- ✓ **No Side Effects**
Prevents burns and ensures safe outcomes


3

Provides ease of use for operators with user-friendly ergonomics.



Tilted tip head moves back and forth, up and down for greater treatment flexibility


- ✓ **Full shot delivery, No missed shots**
- ✓ **Full coverage of the face**
- ✓ **Faster treatment**

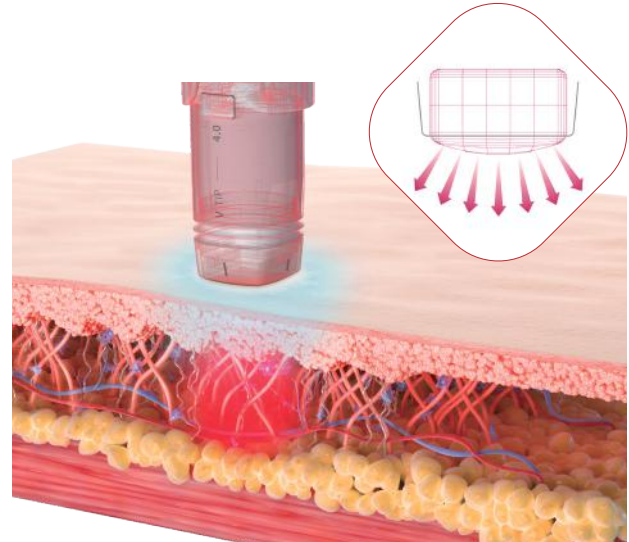

Curved tip surface with radial RF delivery

4

The Hidden Edge Technology

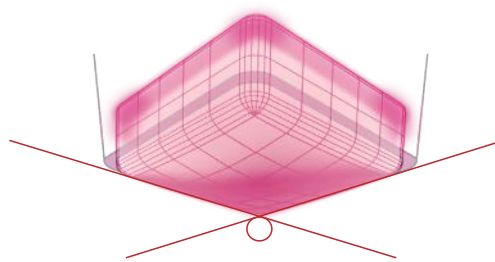

The reengineered tip-edge design eliminates the electrode edge effect, preventing energy concentration at the edges and delivering uniform RF energy across the skin.


- ✓ **Patient safety,**
ensuring no adverse events
- ✓ **Even and consistent energy delivery,**
results in clinically efficacious outcomes
- ✓ **Minimal pain** that requires no anesthesia

5 2 Types of Handpieces and 4 Types of Tips suitable to each indications.


- ✓ **Personalized treatments for both face and body,**
tailored to meet each patient's skin type and treatment goals
- ✓ **No Time Lock,**
No time limitation for tip usage



Various procedural designs are available such as face (eye & forehead), Body and Décolleté

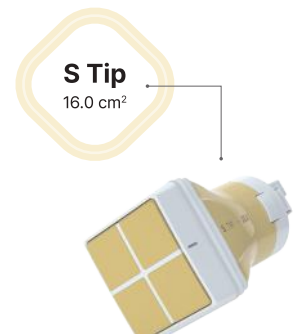

ARMS

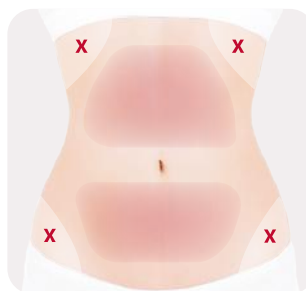

ABDOMEN



KNEES & THIGH




BUTTOCKS

- ✓ One of the few devices with four versatile tips,
enabling tailored treatments without time constraints.



Clinical Efficacy - Improvement in facial skin volume & laxity

Volnewmer effectively treats skin laxity and volume loss with innovative water-cooling technology, ensuring comfortable, safe, and effective results.¹

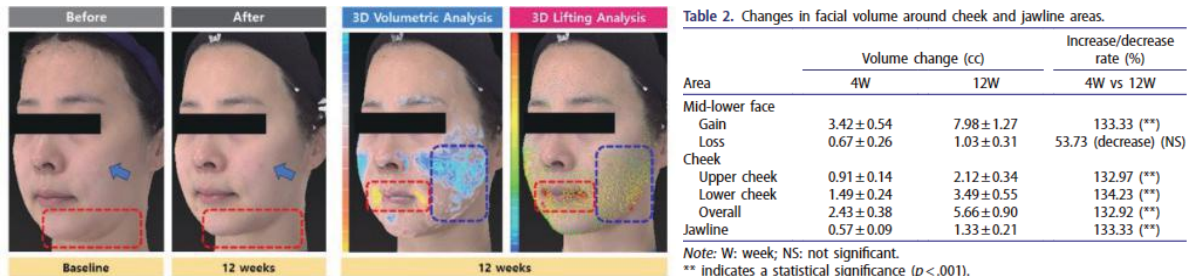

- ✓ The MRF device effectively improves facial skin volume loss and skin laxity, particularly in the cheeks and jawline.
- ✓ The **water-cooling system** ensures moderate pain levels, even without the use of anesthetics, while delivering **high energy** for effective treatment.
- ✓ The study demonstrates that this device is **safe and effective** for treating **mild-to-moderate skin laxity**, especially in **Asian patients**.
- ✓ The **continuous water-cooling system** provides an innovative alternative to conventional methods, addressing their weaknesses by offering better comfort during treatment.

Clinical Efficacy - Reduction of fine wrinkles

Volnewmer's RF therapy reduces fine wrinkles through immediate collagen contraction, utilizing a dual-phase process with optimal energy delivery, while ensuring patient comfort with innovative device tips and techniques.²



- ✓ Volnewmer's RF therapy promotes **immediate collagen contraction** to effectively reduce fine wrinkles.
- ✓ **Dual-phase process** : Immediate and long-term skin improvements through thermal stimulation and collagen remodeling.
- ✓ **Optimal energy levels** and controlled endpoint monitoring ensure safe and effective wrinkle reduction treatments.
- ✓ **Innovative device tips** and sliding technique contribute to patient comfort and enhanced treatment results.



Reference

1. Efficacy and safety of a novel monopolar radiofrequency device with a continuous water-cooling system in patients with age-related facial volume loss. JOURNAL OF DERMATOLOGICAL TREATMENT 2024, VOL. 35, NO. 1, 2333028, <https://doi.org/10.1080/09546634.2024.2333028>
2. Evaluating the Efficacy of Continuous Water-Cooling 115-Watt 6.78-MHz Monopolar RF Therapy for Fine Wrinkle Reduction. Plast Reconstr Surg Glob Open. 2024 Feb 26;12(2):e5623.doi: 10.1097/GOX.0000000000005623. eCollection 2024 Feb.


Before & After




Part	Total Treatment	TIP	LEVEL	Shot
Full Face	1 Treatment	F (3.0cm ²)	3.0	650
		V (4.0cm ²)	0.5	150
		I (0.25cm ²)	1.5	600

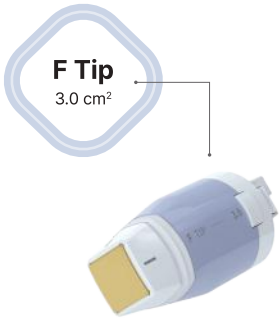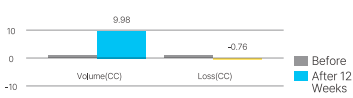


Part	Total Treatment	TIP	LEVEL	Shot
Full Face	1 Treatment	V (4.0cm ²)	2.5	650
		I (0.25cm ²)	2.0	400

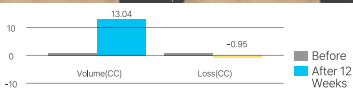

Part	Total Treatment	TIP	LEVEL	Shot
Full Face	1 Treatment	V (4.0cm ²)	4.0	650
		I (0.25cm ²)	2.5	600

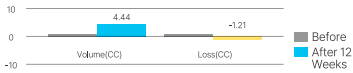

Part	Total Treatment	TIP	LEVEL	Shot
Upper abdomen	2 Treatment	S (16.0cm ²)	6.0	600

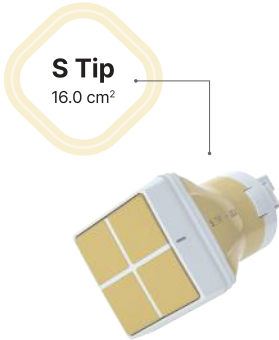

VOLNEWMER



HOME PAGE



Office. +82-2-517-2114 | www.classsys.com | info@classsys.com

All images and texts in this document are intellectual property of Classsys. Any unauthorized usage is strictly prohibited and may result in civil and criminal liability under applicable law. All rights reserved by Classsys.